

Testimonials and Reviews

"We measured the temperature with each of these two thermometers on children and compared the readings to the rectal temperature, which is generally considered the standard. Essentially we found that by every measure that we used, the temporal artery thermometer was more accurate than the tympanic."

- Dr. David Greenes, M.D., Boston Childrens Hospital and Harvard Medical School.

"It is much easier to use than the ear temp probes; no handling of the ear, which typically scares a child, is a great benefit. I have been very impressed with the accuracy of the unit. I would highly recommend it to anyone."

- Dr. Barry W. McCraw, M.D., St. Mary's Regional Medical Center, and Millard Henry Clinic Russellville, AR

"A study published in the Archives of Pediatric and Adolescent Medicine found that the temporal artery thermometer detected fevers more often than the popular tympanic thermometer, which is placed in the ear. The method may even prove to be better than taking a rectal temperature, which has been considered the most practical method of getting an accurate reading."

- Shari Roan, LA Times Health Writer

"The TemporalScanner is the magic wand you've been waiting for."

- Working Mother Magazine

The **EXERGEN**
ComfortSCANNER™
can be purchased at
your local



also available on-line at:

Walgreens.com™

***\$5.00**

Manufacturer's Rebate

Please send this form along with the barcode from the back of the **ComfortScanner™** or **TemporalScanner™** package and original sales receipt to:

**Exergen Corporation
Rebate Offer
400 Pleasant Street
Watertown, MA 02472**

Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Date of Purchase: _____ Where Purchased: _____

*Rebate expires 12/31/08 - max 5 per customer OL1

Questions:

Visit www.temporalscanner.com
Call 1-800-422-3006 toll free

Protected by the following U.S. Patents: 4636091, 5012813, 5199439, 5653238, 5874736, 6045257, 6047205, 6056435, 6292685, 6299347, 6319206, other U.S. and foreign patents pending



EXERGEN Corporation
400 Pleasant Street
Watertown, MA 02472
Tel: 617.923.9900
www.exergen.com



EXERGEN
TemporalScanner™

Recommended for
the entire family!
See Educational Video
www.exergen.com

ComfortSCANNER™
TEMPORAL THERMOMETER



**Accurate
Temperature
with a gentle
Forehead scan**

Method of choice of
premier hospitals

**More
accurate
than ear**
- Harvard Medical
School study
**Accuracy
comparable
to rectal**
- multiple hospital
studies

"Fast, Accurate, Non-Invasive"


- Keith R. Powell, MD - Chairman of Pediatrics
Children's Hospital Medical Center of Akron

**A Better Thermometer
for Your Child**

"I would highly recommend it to anyone."
- Barry W. McCraw, MD

.. Just ask your doctor

Why the Temporal Artery?

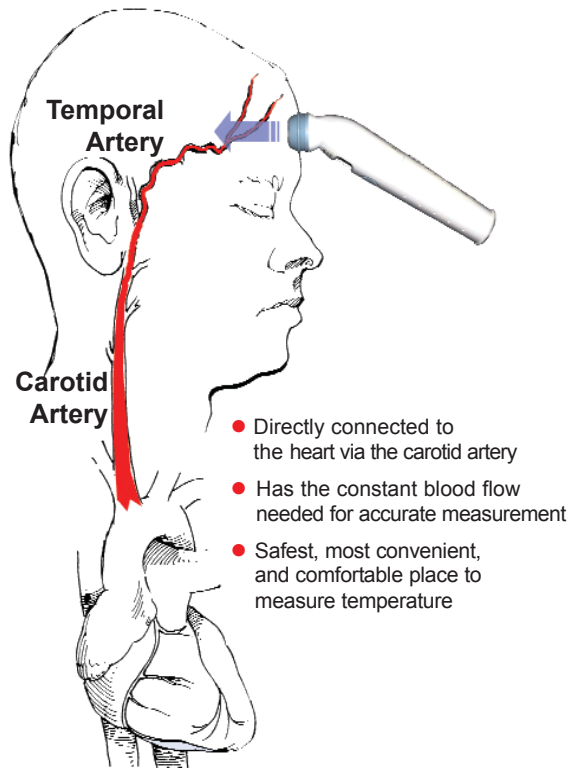

Professional Technology for Parents

...patented **Arterial Heat Balance** provides the breakthrough

- Gentle stroke of forehead captures naturally emitted heat from skin over temporal artery.
- **Red** LED lights and soft beep indicate correct scan
- Takes 1,000 readings per second, selecting the most accurate
- **Arterial Heat Balance** computes the correct temperature



Safe, and Easy to Use

Press and hold button



Lightly scan across forehead

Release button, read temperature

Remove protective cap, read instruction manual completely before using.

The temporal artery is special - the only exposed artery on the skin surface which maintains constant blood flow - necessary for accurate temperature non-invasively at the skin. No other skin area has this property, and no other skin area can give accurate temperatures.

The temporal artery can actually be seen with a high-resolution infrared camera. Note that its location is different for each person, which is why you need to scan laterally to be sure you read it correctly. As the TemporalScanner scans across the artery, the very fast (1000 readings per second) sensing system selects the one peak value, which occurs precisely as the sensor is crossing the artery, and is the most accurate of all. Then the patented Arterial Heat Balance system computes the correct temperature by mathematically replacing the temperature that the blood lost in its short journey from the heart to the forehead.

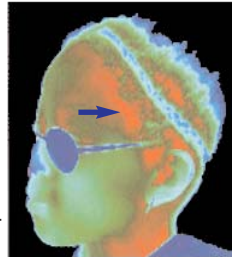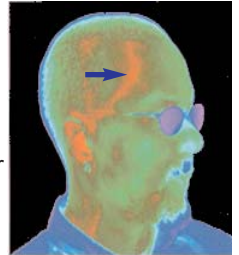

The temperature that you read is called "arterial temperature" and is the most accurate of all temperatures, since it comes from the heart. It will usually be close to rectal temperature, but might be quite different under some conditions, particularly if your child has taken Tylenol, Advil, or other antipyretic, since your TemporalScanner will indicate reduced fever much faster than rectal. Also, when your child develops a fever, the TemporalScanner will indicate it much faster than rectal. As proven in a Harvard Medical School study on infants, your TemporalScanner is more accurate than an ear thermometer; and as proven in multiple hospital studies, has accuracy comparable to a rectal thermometer. This incredible accuracy ... with just a gentle forehead scan.

